


Retinal Gases and Accessories



C₃F₈ Gas

Labtician C₃F₈ gas is used in pneumatic retinopexy procedures and is available in multiuse 450 gram (1.0 lb.) lecture bottles.

Product Code:
MG 66223730

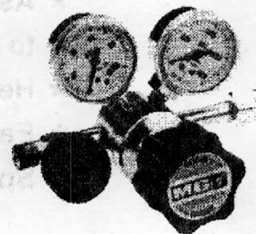

C₃F₈ Regulator

A reusable regulator which allows for the accurate, low pressure delivery of C₃F₈ gas is required. Two gauges are used so that the bottle pressure and the outlet pressure can be read and controlled accurately. (CGA 180 valve)


Product Code:
MG 120-180



SF₆ Gas

Labtician SF₆ gas is used in pneumatic retinopexy procedures and is available in multiuse 227 gram (0.5 lb.) lecture bottles.

Product Code:
MG 3730

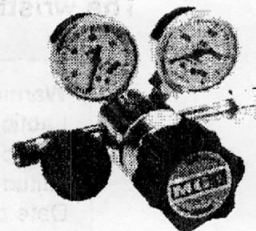

SF₆ Regulator

A reusable regulator which allows for the accurate, low pressure delivery of SF₆ gas is required. Two gauges are used so that the bottle pressure and the outlet pressure can be read and controlled accurately. (CGA 170 valve)

Product Code:
MG 120-170




Accessories

The gas bottle must be upright when used. Labtician offers a wallmount or a table stand for this purpose.

Wallmount
Product Code:
MG 478


Table Stand
Product Code:
MG 477

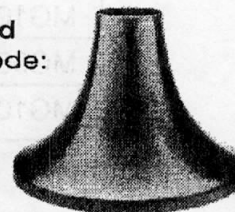

To insure that no gas is wasted, we recommend that the main lecture bottle valve be closed when the bottle is not in use.


Labtician Ophthalmics, Inc.

Unit 6, 2140 Winston Park Drive, Oakville, Ontario L6H 5V5
Tel: 905-829-0055 (Outside of the Toronto area 1-800-265-8391) Fax: 905-829-0056

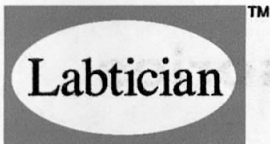

Safety Wristbands For Retinal Gas Patients

- * Assists in explaining to retinal gas patients the need to avoid N2O and high altitude travel
- * Helps ensure patient compliance
- * Each package includes a special marking pen
- * Space for surgeons name and removal date

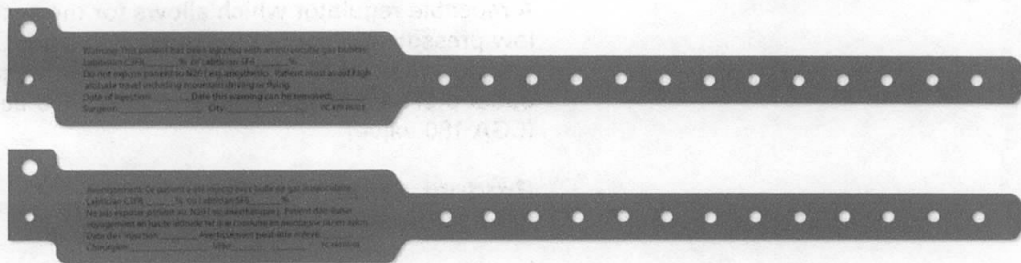

The wristbands are available in both English and French and read:

Warning: This patient has been injected with an intraocular gas bubble:
 Labtician C3F8 _____% or Labtician SF6 _____%
 Do not expose patient to N2O (e.g. anesthetic). Patient must avoid high altitude travel including mountain driving or flying.
 Date of injection: _____ Date this warning can be removed: _____
 Surgeon: _____ City: _____

Avertissement: Ce patient a été injecté avec bulle de gaz intraoculaire:
 Labtician C3F8 _____% ou Labtician SF6 _____%
 Ne pas exposer patient au N2O (ex: anesthésique). Patient doit éviter voyageant en haute altitude tel que conduite en montagne ou en avion.
 Date de l'injection: _____ Avertissement peut-être enlevé: _____
 Chirurgien: _____ Ville: _____

Code	Description
MG500E	Wristbands: 500 pieces with marking pen: English
MG100E	Wristbands: 100 pieces with marking pen: English
MG500F	Wristbands: 500 pieces with marking pen: French
MG100F	Wristbands: 100 pieces with marking pen: French

Labtician Ophthalmics, Inc.
 Unit 6, 2140 Winston Park Drive, Oakville, Ontario L6H 5V5
 tel: 905-829-0055 (outside of the Toronto area 1-800-265-8391) Fax: 905-829-0056